# LIFT mee -Ultrasonic-

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## **INSTRUCTIONS FOR USE**

Congratulations on your purchase.  
We wish you lots of fun with LIFTmee® Ultrasonic!

**IMPORTANT:**

In order to enjoy your LIFTmee® Ultrasonic for a long time,  
read the instructions for use carefully before starting to use the  
device and store it for potential reference later.

**SAFETY**

- Check whether the rated voltage on the type plate matches your local supply voltage. Turn the device off always before you pull out the plug.
- Do not dip device or power pack in water or other liquids.
- Pull off the plug after using the device, before cleaning the unit, or when an error occurs in the power outlet.
- Avoid damage to the power pack by crushing, bending or pulling the cable around sharp edges. Keep the power pack with cable away from hot surfaces and open flames.
- Do not use the device if the power pack or the device itself is damaged.
- Do not leave the device unattended while in use.
- Clean the device carefully with a dry cloth before the first use and after each use.
- Do not use the device for purposes for which it is not intended.
- Do not open the device. No liability is assumed for damages due to improper use.
- Keep the device away from children.
- Use the device only with the supplied power pack.
- This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge, unless they are supervised by a person responsible for their security or they receive instructions from that person on how to use the device.
- Children should be supervised to ensure that they do not play with the device.
- Use the device only with gels and concentrates.

**CONTRAINDICATIONS**

- Cardiac pacemaker
- Pregnancy
- Cancer/Tumours
- Feverish conditions and inflammation
- Bacterial infections of the skin
- Surgical scars within the first 3 months after surgery
- Body areas with severe circulatory disorders
- Circulatory insufficiency
- Coronary artery disease
- Cardiac arrhythmia
- Vascular diseases of the extremities (thrombosis, phlebitis, etc.)
- Coagulation disorders
- Damage to the skin

**BODY REGIONS THAT SHOULD NOT BE TREATED WITH ULTRASONIC/GALVANIC CURRENT::**

- Spinal cord
- Brain
- Eyeballs
- Gonads
- Growth plates in children and adolescents

## LIFTMEE® - ANTI-AGING AT THE HIGHEST LEVEL

We appreciate that you have decided for LIFTmee® Ultrasonic.

With LIFTmee® Ultrasonic, you achieve optimum care results in conjunction with Hyaluron gel and concentrates. Your skin is regenerated from the depths - after a few applications, your skin looks much brighter and fresher. The entirely simple and gentle treatment reduces fine lines and wrinkles, rejuvenates your connective tissue, promotes blood circulation and lymph flow and has a positive effect on skin blemishes and irritations.

## WHAT IS ULTRASONIC/GALVANIC CURRENT?

Ultrasonic is sound waves above the human hearing; they create pressure waves that are used in medicine and physiotherapy. It made a detour and landed in cosmetics: patients ascertained that after ultrasonic testing the treated skin felt smoother and firmer and looked rosier. The effect of ultrasonic on the skin consists of thermal, mechanical and biochemical components. The mechanical effect is achieved by the intensive massage of the elastic fibres of the skin that stimulates the formation of new collagen, by which the outline of the face appears firmer.

In addition, comedones (blackheads) are loosened by the massage and can be removed much easily and painlessly after the ultrasonic treatment. The resulting deep heat in the skin activates metabolism and promotes detoxification, which stimulates the regeneration of the skin and counteracts the aging process. When introducing active agents, they penetrate the skin intensively, so that a visibly more beautiful skin appears in case of rosacea, scars, acne, etc. Particularly high quality Hyaluron gels and concentrates are ideal for the introduction by ultrasound. High and/or low molecular weight hyaluronic acids, which are in high concentration in gel and concentrates, penetrate into deeper skin layers and ensure a lastingly fresh and rejuvenated skin by significantly improving the skin's moisture balance. The purely external use without ultrasound or galvanic current cannot achieve such a remarkable effect.

The ultrasonic treatment is perceived as extremely pleasant; only a slight generation of heat is felt, if at all. You are able to detect a visible improvement of the skin already after 3 to 4 treatments. It also helps to release tension in the neck or relieve headaches. Luigi Galvani published his knowledge in his book "On the Forces of Electricity" in 1791; this knowledge was particularly very well received in naturopathy. Galvanic current over time found its way into cosmetics via the path of natural medicine.

Galvanic current refers to direct current that has no frequencies, i.e. the current flows always in one direction. Positively charged particles move towards the - pole, negatively charged particles migrate to the + pole and induce appropriate responses there; the effect of the treatment is based thereon.

Positive features in the cosmetic field: effective cleaning of the pores, active agents are introduced into the skin more effectively, the connective tissue is strengthened, blood circulation and lymph flow are lastingly improved.

## GIVE YOUR SKIN THE ULTIMATE CARE EXPERIENCE.

The outstanding success of LIFTmee Ultrasonic is based on the perfect combination of very high quality ingredients and innovative technology. LIFTmee Ultrasonic is responsible for transporting the ingredients. Hyaluron gel and the concentrates supply the highest quality ingredients. A perfectly coordinated team that optimally provides your skin a fresher and more radiant look in the long term. You will notice visible results after only a few applications. Select the individual treatment programme that is good for your skin. Clean your face as usual before each treatment.

By the way, you not only please your facial skin with an increase in care, but also make your décolleté or hands look fresh without fuss. And even with cellulite, an ultrasonic treatment on the affected parts of the body can achieve obvious results.

## WHAT EXACTLY DO YOU EXPECT?


After the first treatment, your skin is better supplied with blood; it looks radiant and rosier with increased supplies of oxygen. At the same time, waste products are removed effectively, which you can detect, among others, in the reduced under-eye bags. Thanks to the deep massage your connective tissue is trained, creases and fine lines disappear and the depth of the wrinkles significantly softened. Your face in general looks brighter and fresher with defined contours. Likewise, the regeneration of cells is activated effectively - your skin gets more tone.

Thanks to the disinfecting effect of ultrasonic, skin blemishes are significantly reduced.

## CHARGING

The device has to be fully charged before it is used the first time. The standard recharging time is approx four hours. When the device is fully charged, the display light turns green.

## FIXING THE COTTON PAD FOR CLEANSING
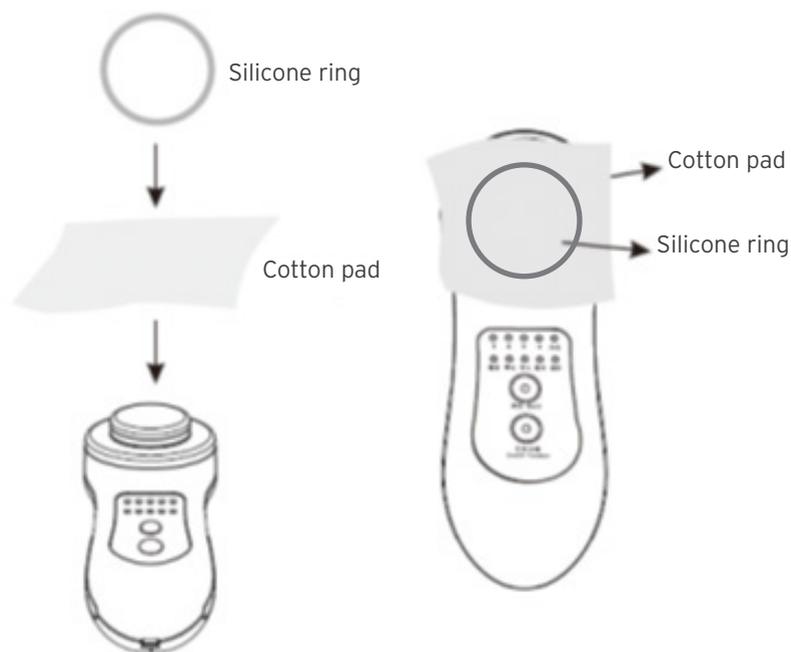
## CLEANING PROGRAMME

Remove your make-up.

Switch on (ON/OFF) the device and select the programme „clean“ also select the desired intensity by pressing the “+” button.

Generously apply a cleansing tonic or just water on the ultrasound head or on the skin to be treated and place the ultrasound head on the skin area to be treated. Guide LIFTmee - Ultrasonic - with slow, gentle, circular movements - without exerting pressure - over the skin. Make sure that the ultrasound head is rotating always on a well moisturised skin surface.

If necessary, apply cleansing tonic or water again.

After 10 minutes the device switches off automatically, and you will be alerted by audible beeps. In this programme, LIFTmee - Ultrasonic - emits positive ions, which means negative ions are extracted from the surface tissue of the skin, cleaning the deep skin effectively. Depending on how your skin is inclined to blemishes, treat yourself with the cleaning feature regularly or sporadically.

## PENETRATION OF ACTIVE INGREDIENTS

Switch on (ON/OFF) the device and select the programme „IN“.

Also select the desired intensity by pressing the “+” button.

Generously apply Hyaluron gel or Hyaluron concentrates on the ultrasound head or on the skin to be treated and place the ultrasonic head on the treatment area. Guide LIFTmee - Ultrasonic - with slow, gentle, circular movements - without exerting pressure - over the skin. Make sure that the ultrasound head is rotating always on a well moisturized skin surface.

If necessary, apply serum or gel again.

After 10 minutes the device switches off automatically, and you will be alerted by audible beeps. In this programme, the active agents of hyaluron gel or hyaluron concentrates are introduced deeper into the skin due to the migration of ions. You should treat yourself with this luxury care in combination with the ultrasound deep massage 5 twice a week regularly to draw the utmost benefit from the valuable and highly effective ingredients.

## ULTRASONIC DEEP MASSAGE

Switch on (ON/OFF) the device now you are in the Ultrasonic Modus.

Select the desired intensity by pressing the “+” button.

Generously apply hyaluron gel or concentrates on the ultrasound head or on the skin to be treated and place the ultrasound head on the skin area to be treated.

Guide LIFTmee - Ultrasonic - with slow, gentle, circular movements - without exerting pressure - over the skin. Make sure that the ultrasound head is rotating always on a well moisturised skin surface.

If necessary, apply serum or gel again.

After 10 minutes the device switches off automatically, and you will be alerted by audible beeps. In this programme, the active ingredients of hyaluron gel or concentrates are introduced deeper into the skin by the deep massage, simultaneously causing a kind of training effect on the connective tissue that allows for better blood circulation and the effective removal of waste products via the lymphatic system. In addition, the pores of your skin will be cleaned due to the disinfecting effect of the ultrasound.

You should treat yourself with this luxury care in combination with the penetration of active agent ION - twice a week regularly to draw the utmost benefit from the valuable and highly effective ingredients.

## LIGHT THERAPY

Clean the face.

Switch on (ON/OFF) the device and select the programme “Light”.

Select the desired light by pressing the “+” button.

### MODE 1 = RED LIGHT (630NM)

Used as a long-wave red light, collagenlight is helpful for the skin's renewal process.

The colour of blood and heat penetrates deeply into the skin until the layers, where the collagen and elastin can be found. Collagen and elastin are structure proteins that provide our skin with firmness and elasticity. Collagen light supports the elimination of the hardenings and deposits of these proteins.

It additionally strengthens the skin's stability and elasticity.

### MODE 2 = ORANGENE LIGHT ( 610NM)

Activates the cell's microcirculation and provides an even and bright complexion.

### MODE 3 = BLUE LIGHT (463NM)

The colour of water and sky calms irritated skin and has an antiseptic effect.

For this reason, blue light is particularly suitable for improving the skin appearance in case of blemished skin and couperose. Additionally, it stimulates the skin's oxygen production, cleanses and combats the new formation of bacteria. The skin's immune system can be improved.

Generously apply hyaluron gel or hyaluron concentrates on the ultrasound head or on the skin to be treated and place the ultrasound head on the treatment area.

Guide LIFTmee Ultrasonic with slow, gentle, circular movements and without exerting pressure, over the skin.

Make sure that the ultrasound head is rotating always on a well moisturized skin surface. If necessary, apply serum or gel again.

## VIBRATION MESSAGE

We recommend to use this mode in a combination with following functions:

Cleansing, Penetration of active ingredients and ultrasonic.

The soft vibration stimulates the lymph, reduces dark circles under the eyes and tightens the facial contour.

Clean the face. Switch on (ON/OFF) the device and select the "Vibration" programme by pressing the button repeatedly.

Also select the desired intensity by pressing the "+" button.

Generously apply hyaluron gel or hyaluron concentrates on the ultrasound head or on the skin area to be treated and place the ultrasound head on the treatment area. Guide LIFTmee - Ultrasonic- with slow, gentle, circular movements and without exerting the pressure, over the skin. Make sure that the ultrasound head is rotating always on a well moisturized skin surface. If necessary, apply serum or gel again.

## COMBINED FUNCTIONS

In addition to the individual functions the LIFTmee Ultrasonic offers you following combination possibilities:

- 1) Ultrasonic + Cleansing + Vibration
- 2) Ultrasonic + Penetration + Vibration
- 3) Ultrasonic + Vibration
- 4) Ultrasonic + Cleansing + Light therapy + Vibration

## ENVIRONMENTAL INFORMATION

This product is made of high quality parts and materials.

At the end of its service life, do not dispose of the product with the normal household waste, take it to a collection point for recycling electronic and electrical equipment.

Help to protect the environment by recycling used products.

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Power supply: In AC110-240V, out DC 5V

Output: 1.11 W

Weight: 150g

Ultrasonic: 3 MHz

Technical and optical modifications reserved.

Do not disassemble the device under any circumstances.

No liability is assumed for damages due to improper use.

Do the function test - you can see the ultrasound vibrations for yourself by the use of a drop of water. Moisten the ultrasound head with a drop of water, and switch in the Ultrasonicmode.

You will see how the water drop vibrates on the ultrasound head. After the treatment - depending on skin type - you may experience a feeling of tightness; remove residual product from the skin like a peel-off mask. Therefore clean your face with a damp cloth after every treatment and then apply hyaluron gel or hyaluron concentrates sparingly on the facial skin. Continue with your usual day or night cream.

Optimal results are achieved when you use hyaluron gel or hyaluron concentrates on a regular basis as a foundation under your normal day and night cream. In addition treat yourself twice weekly with LIFTmee®, and after few applications you will be looking much more relaxed in your daily life.

You would not want to miss LIFTmee® any more!

**WM BEAUTYSYSTEMS AG & Co. KG**

Carl-Friedrich-Gauß-Straße 11

50259 Pulheim

Tel +49 (0)2234 - 98422 - 0

Fax +49 (0)2234 - 98422 - 711

Internet [www.wm-beautysystems.com](http://www.wm-beautysystems.com)

E-Mail [info@wm-beautysystems.com](mailto:info@wm-beautysystems.com)